



Minimum height: 19cm.




Maximum height: 68cm.

Improves patient position and enhances their recovery

The cardiac chair position is essential for patients with heart and breathing problems, helping to enhance their recovery.

Position ergonomy: auto-regression

The regression function of MEDI LOW bed allows movement in the backrest (11cm), reducing friction to the minimum.

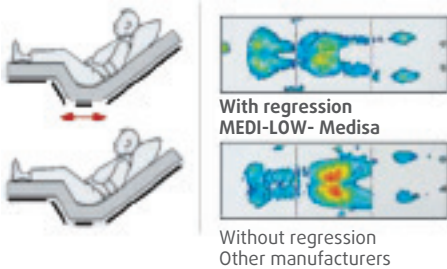

MEDI-LOW Bed

Comfort and individual care thanks to its multiple functions

Basic Equipment

- 1. Fully electric bed with scissor-lift raising system.
- 2. Maximum safe working load: 230 Kg.
- 3. Bed mattress made up of four independent ABS modules.
- 4. Patient handset with function cancellation.
- 5. Backrest regression (11 cm).
- 6. 75 mm diameter wheels.
- 7. Supports for IV pole and lifting pole.
- 8. Mattress retainers.
- 9. Central braking system.
- 10. Bilateral accesories holder.
- 11. 8 hooks for restraining straps.

Accessories

- 1. Urine bag holder (REF. CHA20).
- 2. Lifting pole (REF. CHA17).
- 3. IV pole (REF CHA15).

Optional equipment

- 1. Steel side railles (REF. CHV50).
- 2. Lateral side rails in vaporized beech (REF. CHV56).
- 3. Headboard/footboard made in vaporized beech(REF. CHV56BG/CG).
- 4. Melamine headboard/ footboard (REF. CHV56B/C).
- 5. Emergency battery (REF. 566).

Electric features

Power supply	230-110V /50-60 Hz
Maximum current intake	1.6A, 370 Volt
Protection indicator	IP 66
Class of protection	Class II
Electric shock protection	Type B

Technical features

Overall dimensions	206 x 104 cm
Patient surface dimension	200 x 86 cm
Height adjustment	19- 68 cm
Backrest angle	70°
Legrest angle	30°
Backrest regression	11 cm
Trendelemburg/reverse trendelemburg	±16°
Weight without accessories	150 kg
Maximum safe working load	230 kg



www.medical-iberica.com

Lanzahita, 6 - Pol. Ind. Albarreja - 28946 Fuenlabrada, Madrid (Spain) - Ph. +34 91 208 94 00 - F. +34 91 208 94



Hospital furniture
MEDI-LOW Bed

Comfort and individual care
thanks to its multiple functions

REV. 12/2017 NOTE: Measures and features of the presented products may be modified without prior notice.

MEDI-LOW Bed

Comfort and individual care thanks to its multiple functions



Reduces physical effort needed to move patients.

Helps to the recovery of pressure bedsores.

Innovative concept of lateral side rails which cover the 100% of the mattress, according to EN 60601-2-52 rule.

With the new systems incorporated into the MEDI-LOW bed, we help to speed up the patient's recovery thanks to an advanced ergonomic design that increases their autonomy and mobility.

New safety requirements

According to EN 60 601-2-52 standard for electric beds, side rails should minimize the risk of trapping the patient's head, neck or limbs in order to avoid possible fatal accidents.

Steel side rail option



Positioning and comfort



Patient handset with functions cancellation

Easy cleaning

Our high quality products, reduce time for cleaning up to 30%. Smooth surfaces and detachable panels allow a quick and easy cleaning and disinfection.



Trendelenburg $\pm 16^\circ$



Reverse trendelenburg $\pm 16^\circ$

Ergonomic head and footboard panels that help manoeuvring the bed. The materials used (varnished beech) meet all the necessary requirements for hygiene and cleaning.

TOTAL PROTECT side rails: They offer perfect protection against falls. Compliant with EN 60 601-2-52. Made in vaporised beech.

Four easily detachable ABS sections.

Bilateral accessories holder.

Supports for fixing straps.

Headboard and footboard made in vaporized beech.

Central brake. Rubber wheels with retractable system